

memmert

trust the best



Waterbath

WTB

accessories

WTB | PRACTICAL AND ROBUST

Easy to clean

- rounded corners
- stainless steel surfaces

Easy to empty

- new drain system

Easy to handle

- smallest possible footprint




WTB | INTUITIVE AND USER-FRIENDLY

Easy to operate

- 3.5 inch touchscreen
- New user interface

Versatile functions

- Timer
- Alarm
- 2-point calibration
- Setting Celsius / Fahrenheit

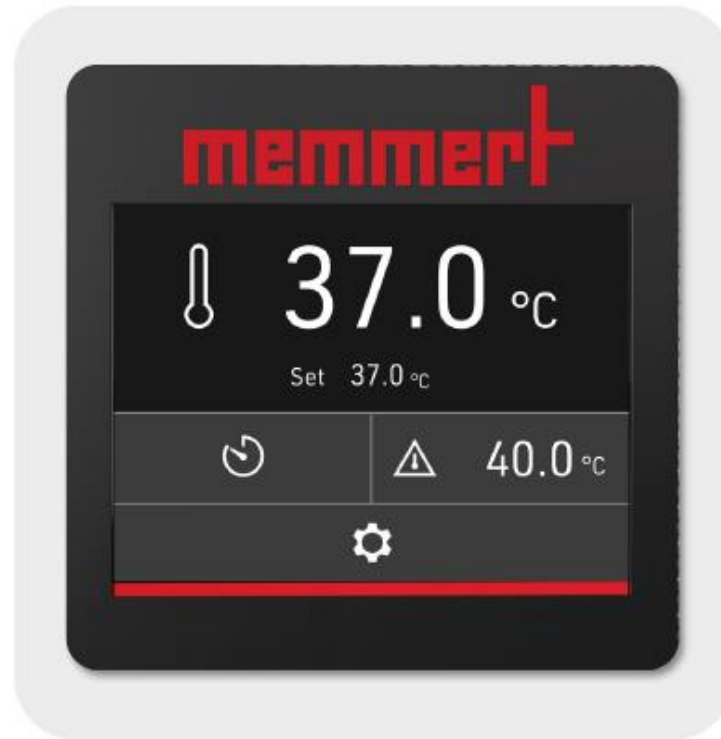

WTB | CUSTOMVIEW

Thanks to the customisable remote View CustomView you always keep an eye on the important information.



WTB | TECHNICAL DATA



WTB6

7,5 l

WTB11

10 l

WTB15

15 l

WTB24

24 l

WTB35

37,5 l





WTB50

51 l

Model sizes / Article description	6	11	15	24	35	50
Power consumption at 230/115 V, 50/60 Hz (W)		1000		1420		2800
Operating temperature range (°C)	min. 5 above ambient temperature up to +100					
Setting temperature range (°C)	+10 to +100					

WTB ACCESSORIES | OVERVIEW

Flexible and modularly expandable - adapted to your needs.

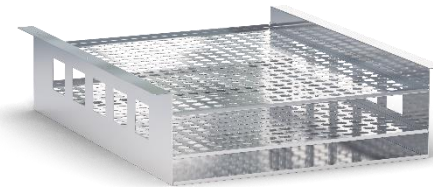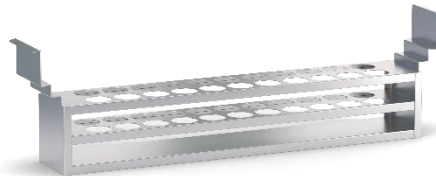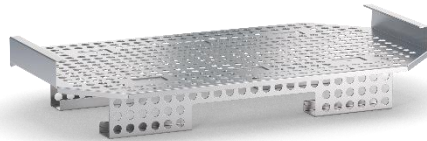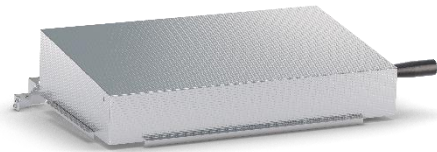

WTB ACCESSORIES | SHAKING DEVICE SV

Maintains load movement.

- Two sizes: SV1524 and SV3550.
- A shaking frequency 30 - 150 strokes per minute to freely choose from
- Adjustable program runtime and alarm protection
- Storage basket and insert racks are present to prevent the load from sliding or tipping over



WTB ACCESSORIES | PELTIER-COOLING DEVICE CDP115

Compact design and cooling with environmentally friendly Peltier technology.

- One cooling device for all sizes: CDP115 and CDP115X2
- Cooling capacity up to 10 °C with Peltier technology
- Both environmentally friendly and energy efficient



Model sizes / Article description	6	11	15	24	35	50
Cooling rate (K/h)	14,1	9,9	7,1	4,5	3,4	2,2

WTB ACCESSORIES | PUMP DEVICE CPM650



Improvement of the thermal homogeneity.

- A continuous circulation of the temperature control medium
- Self-priming pump is suitable for continuous operation
- Flow rate of 1.2 liters per minute



WTB ACCESSORIES | REVERSIBLE SHELF AND DIFFERENT COVER VARIANTS

- Protection against foreign bodies
- More efficient and energy-saving operation
- Reduction of evaporation
- Uniform heating of the temperature control medium

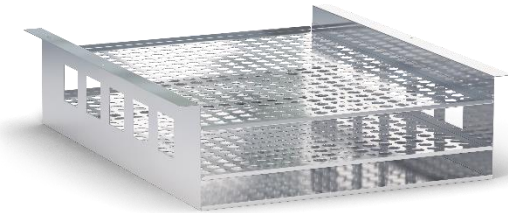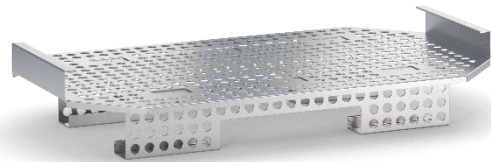

WTB ACCESSORIES | APPLICATIONS

- Temperature control of ointments, emulsions, samples and nutrient solutions
- Sample preparation for photometric examinations
- Investigations of soil samples
- Enzymatic rapid tests for the detection of heating in milk
- Beer enforcement tests
- Alternating tests for resistance testing of O-rings
- Heating of cell media in the laboratory
- Warm storage for material relaxation of components
- Calibration of articles (e.g. bath thermometers) at different water temperatures
- Sous-vide cooking meat and fish
- ...